
DMC-EZ (DICOM Media Creators)

Web/CD Viewer User's Guide

For DMC-EZ Series 2.x Products

Ver. 2.2

October 2012



SST Group Inc.
309 Laurelwood Rd. #20
Santa Clara, CA 95054
800-944-6281
408-350-3450
408-350-3100 (fax)
www.sstgroup-inc.com

Table of Contents

DMC-EZ (DICOM Media Creators)	1
Section 1 Viewer (CD/WEB).....	3
1.1 Installation of EzDicom Web Viewer	3
1.2 Quick Start.....	4
1.3 Common Tools.....	4
1.3.1 Pan.....	4
1.3.2 Windowing.....	4
1.3.3 Annotations and Measurements	5
1.3.4 Zoom	5
1.3.5 Original	5
1.3.6 Negative	5
1.3.7 Orientation.....	6
1.4 Layouts	6
1.4.1 Series Layout	6
1.4.2 Images Layout.....	7
1.5 Other Tools.....	7
1.5.1 Comparative Mode.....	7
1.5.2 Groups.....	8
1.5.3 Video Control.....	8
1.5.4 Ref Lines.....	8
1.5.5 Magnifier	8
1.5.6 Image Infos.....	9
1.5.7 Current / All Series	9
1.5.8 DICOM Header.....	10
1.6 Export	11
1.7 Print	12
1.8 Report.....	13
1.9 Key Images.....	13
1.10 MPR	13

Section 1 Viewer (CD/WEB)

1.1 Installation of EzDicom Web Viewer

To install EzDicom Web Viewer on a remote computer, open your internet browser and navigate to DMC-EZ login page. (IP address of DMC-EZ)




A blue rectangular login form titled "Connection to DMC-EZ". It contains two input fields: "User Name" and "Password". Below these fields is a "Login" button. In the bottom right corner, it says "Version 2.0.2".




Click on the link on the right and select Run to start the installation. Follow the instructions provided in the installation wizard.



1.2 Quick Start

Insert the CD in the drive. If Autorun is enabled, EzDicom CD Viewer will start automatically. If not, browse the CD folder and double click on Launcher.exe. Or if you are connected to the DMC-EZ you can launch the Ez-Dicom Web Viewer by selecting studies.  and click "Visualize" or double click study.

The first Study will be loaded when EzDicom CD Viewer starts.
You can navigate through Studies and Series using the Series panel on the left:
You can also navigate using the navigation buttons in the toolbar:

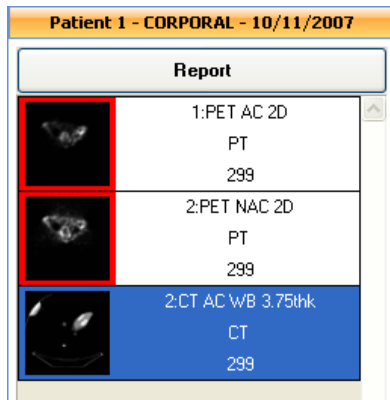

1.3 Common Tools


1.3.1 Pan

The Pan tool allows you to move the Image in the viewer. Left click and drag on the Image to Pan.



1.3.2 Windowing

The Windowing tool allows you to modify the Level and Width of the current Series. Left click and drag on the Image to modify those values. You can also use a Predefined Windowing. Click on the arrow on the right of the Windowing button and select a Preset.

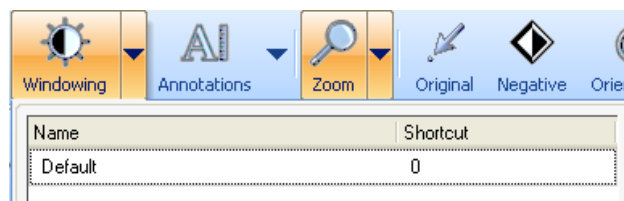

1.3.3 Annotations and Measurements


The annotations and measurements tools allow you to quickly add measures and annotations to Images. Click on the Annotations button and select the measure or annotation you want to add.



1.3.4 Zoom

You can zoom in the current Series using one of the Zoom tools.

With the Normal Zoom tool, you can simply zoom in or out the Image with a right click and drag. The Pan and Zoom tool will also center the Image where you clicked.

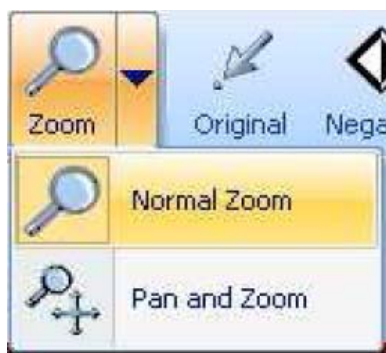

1.3.5 Original

The Original tool returns the current Series to its original windowing, zoom, rotation and position.



1.3.6 Negative

The Negative tool allows you to invert the colors for the current Series.



1.3.7 Orientation

You can apply rotations or flips to the current Series using the Rotation tools. You can rotate an Image clockwise or counterclockwise, as well as flip the Image horizontally or vertically.



1.4 Layouts

1.4.1 Series Layout


The Series Layout tool allows you to visualize multiple Series at the same time.



Click on the Series button and select a Layout. You can then add Series on each layout by selecting the layout area on the viewer and clicking on a Series in the Series panel.

1.4.2 Images Layout

The Images Layout tool allows you to visualize multiple Images from the same Series at the same time.

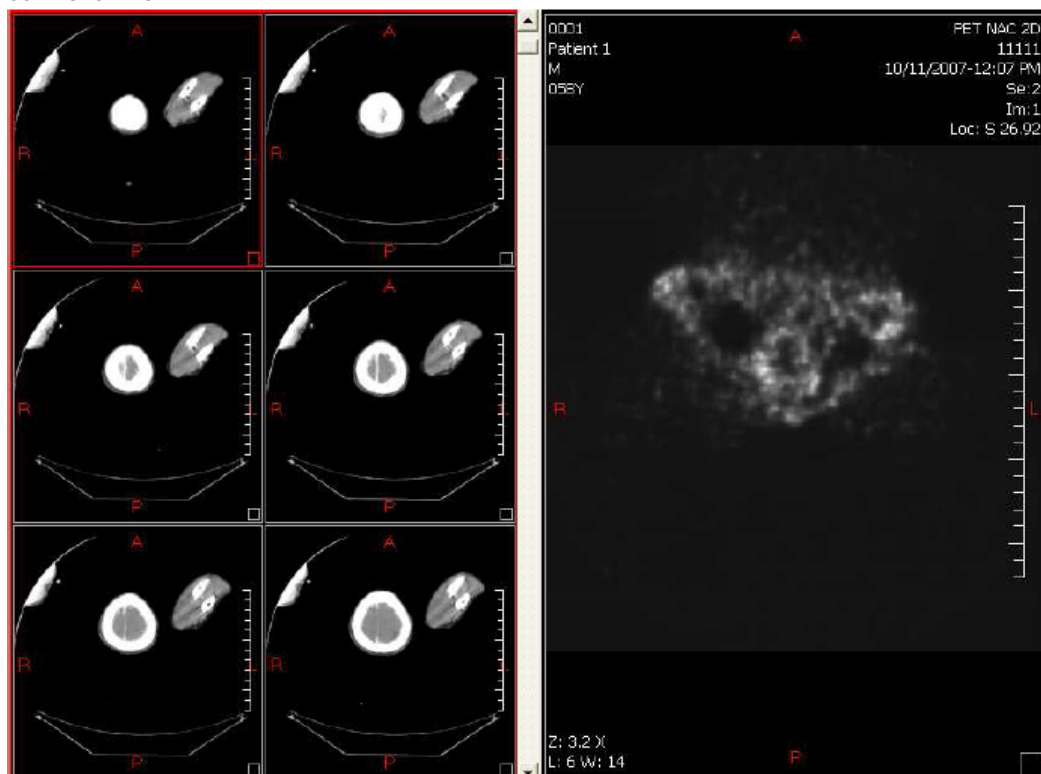

Simply click on the Images button and select a Layout.

1.5 Other Tools


1.5.1 Comparative Mode

The Comparative Mode allows you to scroll into different Series at the same time. Select a 2x1 Series Layout, add a Series in each Viewport, Click on comparative Mode and scroll on a Series. Both Series will scroll at the same time.



1.5.2 Groups


The Groups tool is an extension of the Comparative Mode. You can assign Series to a group. If you do so, all the Series of this group will be in Comparative Mode.

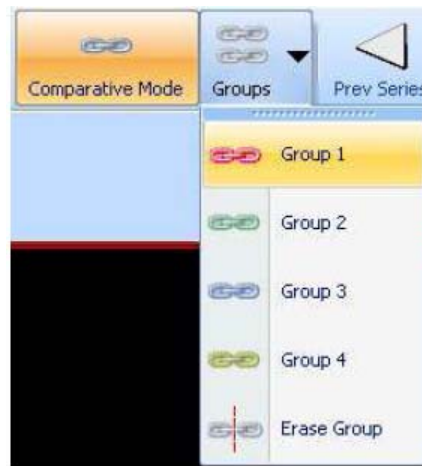

For example, assign two or more Series to a group and scroll. All the Series assigned to this group will scroll together. You can create up to 4 groups.

1.5.3 Video Control

With the Video control tools you can visualize Series like a film, and you can choose the Cine speed and its direction.

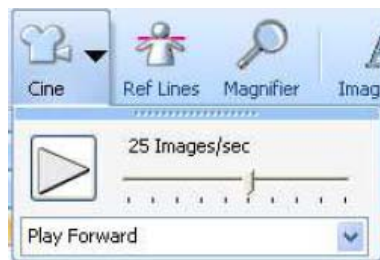

Click on the Cine button to access Cine tools.

1.5.4 Ref Lines

The Reference Lines tool shows the location of the active Image as Lines in the other Series of the Study. Simply enable this functionality by clicking on the Ref Lines button.



1.5.5 Magnifier

A Magnifier is available by clicking in the Magnifier button.



You can move the Magnifier window by left clicking and dragging it. You can also increase the magnifier window size and zoom.


1.5.6 Image Infos

The Image Infos tool allows you to toggle information display on the Image.



1.5.7 Current / All Series


The Current Series / All Series tool defines whether the tools will modify the current Series or all the Series displayed in the viewer.

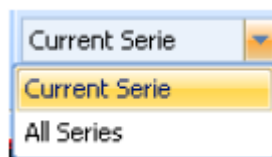

For instance, select All Series in the list modify the windowing of a Series and all the displayed Series will be affected.


1.5.8 DICOM Header

Using the DICOM Header you can easily access the DICOM information contained in the Image. Right click on the Image and select the DICOM Header tool.

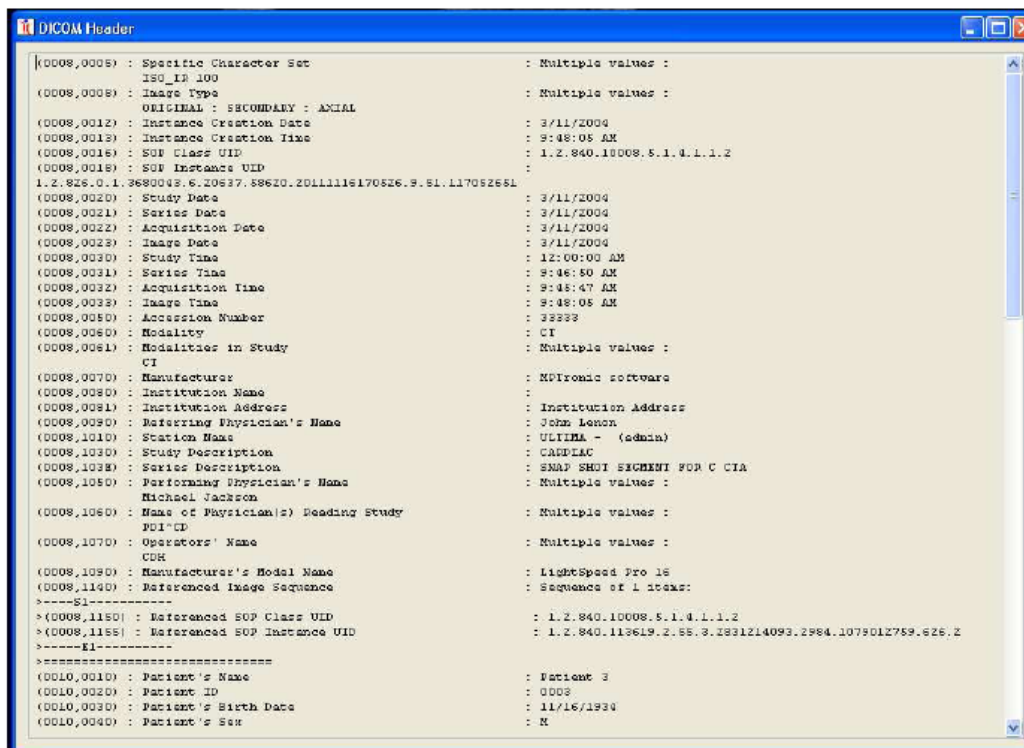

1.6 Export


With the Export functionality you can export Images in several formats. Click on the Export button to open the Export window.

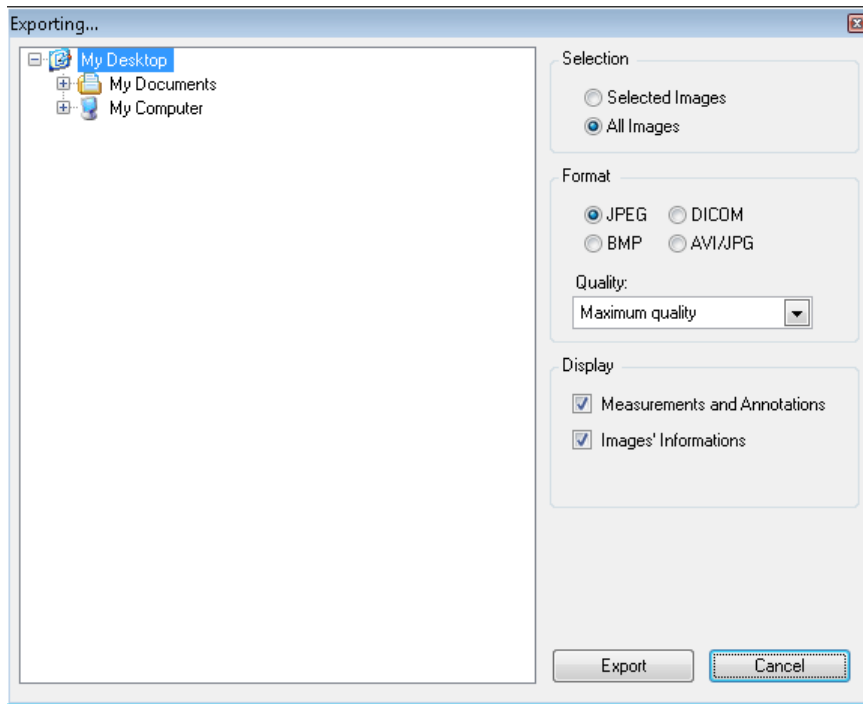

On this window you can choose to export Images with or without Image Information, measurements and annotations, as well as choosing the exportation format and quality. You can also choose whether you want to export only the current Series or all the current Study.

1.7 Print


The Print tool allows you to print selections of Images in multiple layouts.

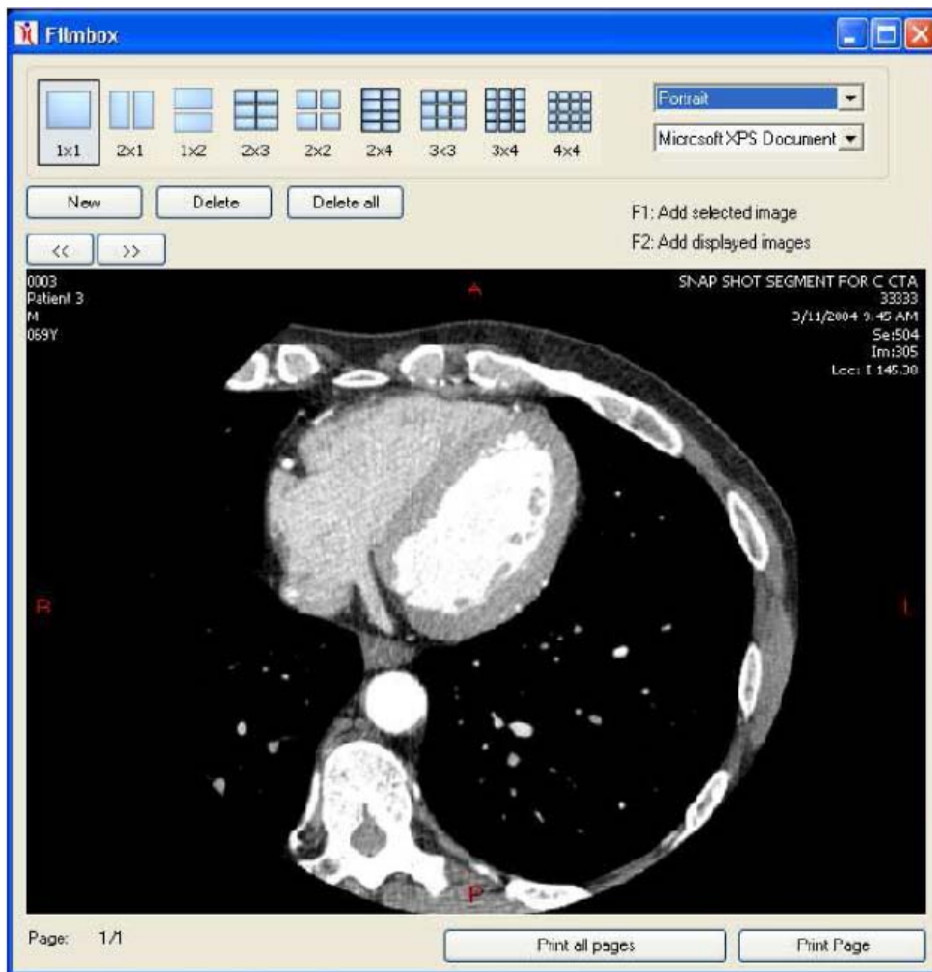

You can add selected Image to the filmbox by pressing F1, or all the displayed Images by pressing F2. Once you have added all Images you want to the filmbox, click on Print to open the Print window. You can there select the layout you want and print Images.

1.8 Report


You can edit or add a Report to a Study using the Report window. Each Study loaded in the Viewer has a Report button, click this button to open the Report window.

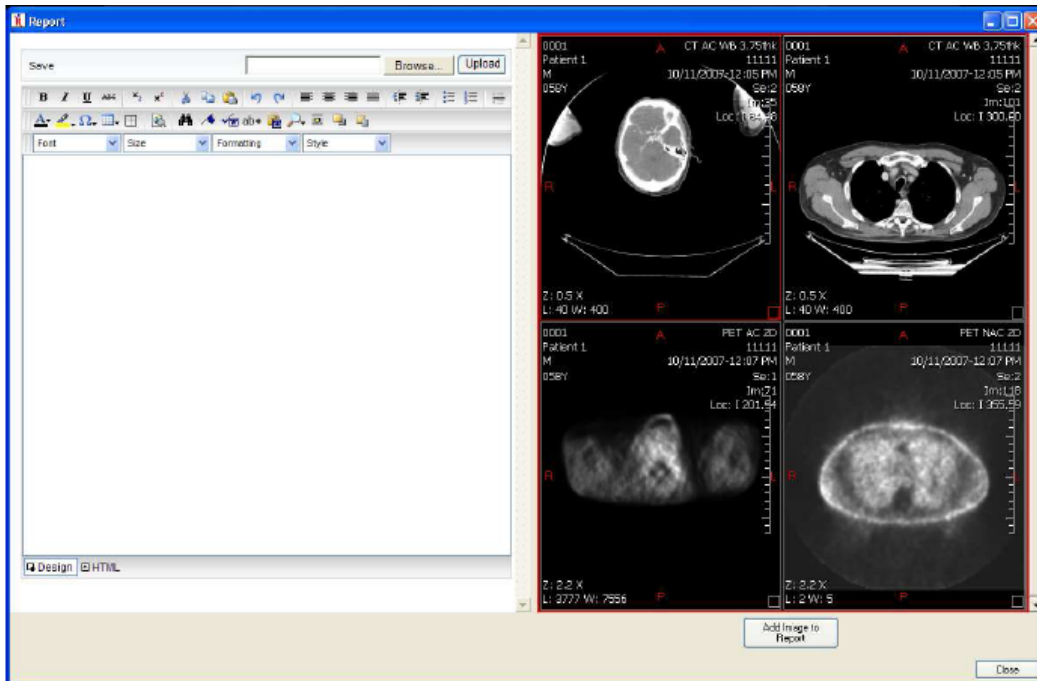

In this window you can type in the report using the basic text processing tools. You can also add links to some Images to the report. Select an Image on the Viewer on the right of the window and click on Add Image to Report. To have Images available in the Report window, you need to add them first by selecting them in the viewer and pressing F3. Once you are done typing the Report, click on Save, and then on Validate to validate the Report.

1.9 Key Images

You can easily visualize a Study Key Images using the Key Images window. If a Study has Key Images available, click on the Key Images button.

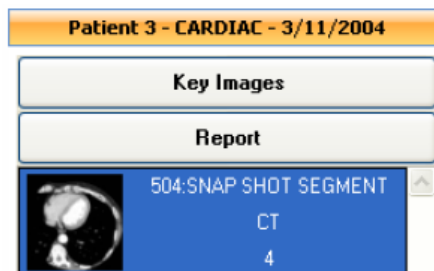

1.10 MPR

EzDicom CD Viewer MPR is a separate module and has its own documentation.